
FeelGood2

(Series JA)

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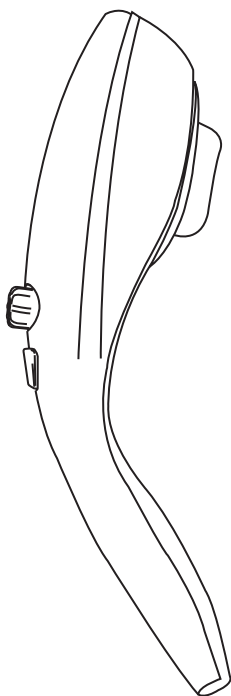Адрес производителя: D Building No. 100, JINFU ROAD, Chengnan Industrial Zone, TONG'AN District, Xiamen, China

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

handheld massager

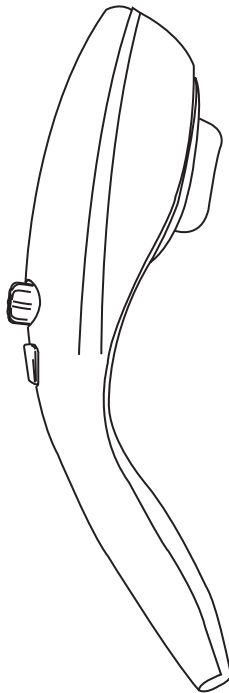User Manual

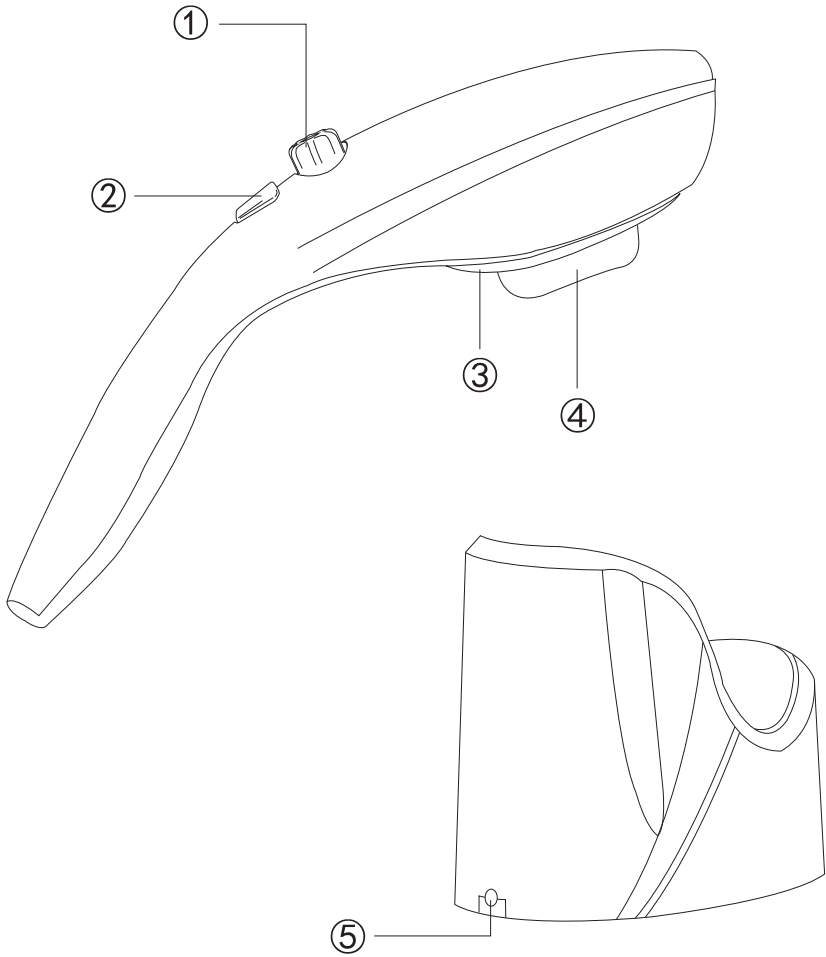
- ① Power / Speed Knob with Indicator
- ② Infrared Button
- ③ Infrared Light
- ④ Vibration Massager
- ⑤ Power Port


Instructions

1. Power On, Power Off

Turn the Power / Speed Knob clockwise / counterclockwise with a click sound to turn on / off the massager, the indicator will turn on / off accordingly. The default massage cycle is 15 minutes after which the massager will turn off.

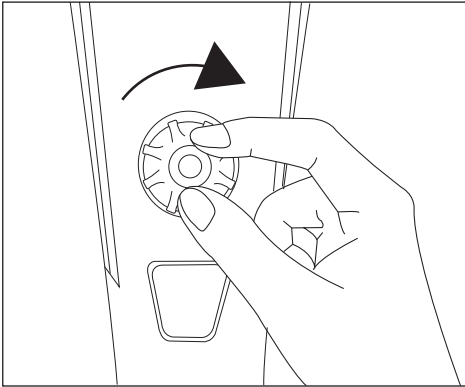
Power On

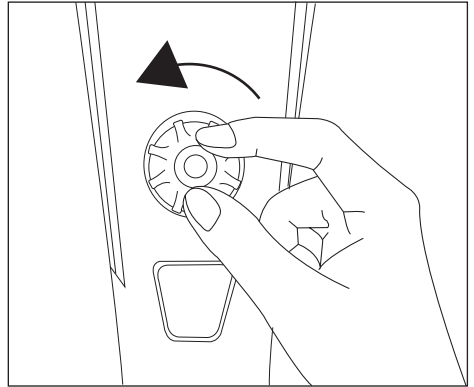
Power Off

2. Speed Adjustment

When the massager is turned on, turn the Power / Speed Knob clockwise / counterclockwise to increase / decrease the messaging speed.


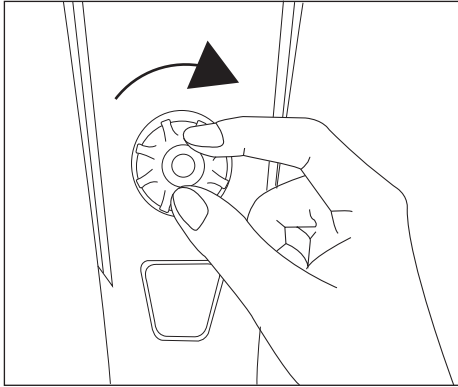
High Speed

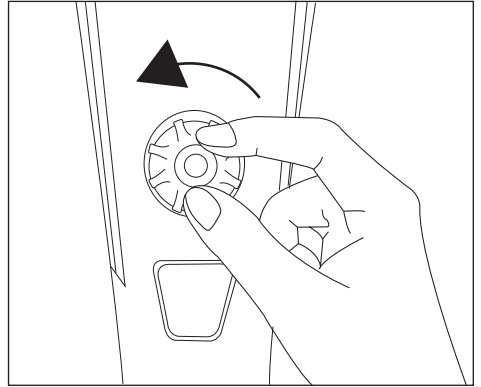Low Speed

3. Infrared Heating

When the massager is turned on, press the Infrared Button to turn on / off infrared heat, the infrared light will turn on / off accordingly. The infrared massage function can also be turned off with the power.


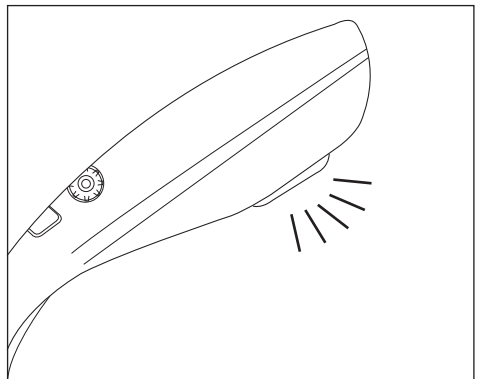
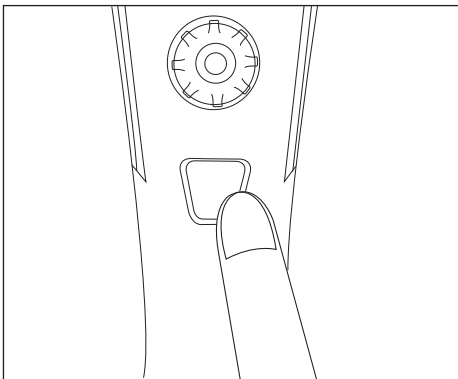
Charging the Massager

1. Place the massager into the charging dock.
2. Plug the power cable into the power port.
3. Connect the adapter to a power source.
4. The power indicator lights up blue when charging is started.
5. The power indicator turns automatically when fully charged after 3.5 hours.

Note:

- When the battery is low and the massager is turned on, the power indicator will flash.
- Do not turn on the massager when it's being charged in the dock.

Massage Areas

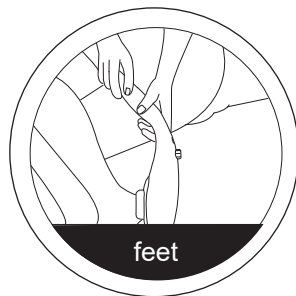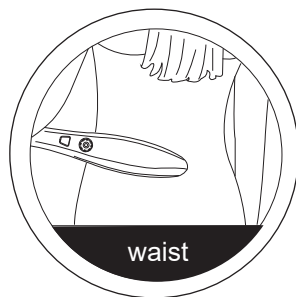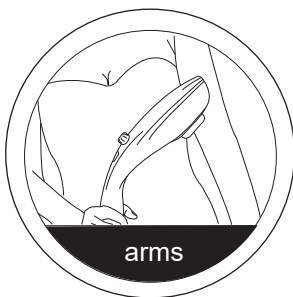Important Safety Instruction

When using the product, all instructions and basic precautions should always be followed, including the following:

DANGER - To reduce the risk of electric shock: Always unplug this product from the electrical outlet immediately after use and before cleaning

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) The product should never be left unattended when plugged in. Unplug from outlet when it is not use and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Never operate the product with the air openings blocked. Keep the air openings away from lint, hair, and something like that.
- 4) Do not drop or insert any object into any opening.
- 5) Do not use outdoors.
- 6) Do not operate where aerosol (spray) products are being used or where oxygen is administered.
- 7) Be sure to pull the plug out of the socket when you reach in. And don't pull it out when your hands are wet. That could cause electric shock and injury.
- 8) Do not carry this product by its supply cord or use cord as a handle.
- 9) Keep the cord away from heated surfaces. If the fault is found, stop using immediately and cut off the power supply. Avoid accidents or injuries. please do not repair it by yourself, please contact the seller or manufacturer in time, avoid accidents or injuries.
- 10) Remove the power plug from the socket immediately after the power failure. Otherwise, it may cause electric shock, short circuit and fire.
- 11) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examinatio.
- 12) Special attention to components most susceptible to wear an advice to replace defective components immediately and/or keep the equipment out of use until repair.
- 13) Use this appliance only for its intended use as described in this manual.
- 14) The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- 15) The appliance is only to be used with the power supply unit provided with the appliance.

16) This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.

※ Children should be supervised to ensure that they do not play with the appliance.

※ Cleaning and user maintenance shall not be made by children without supervision.

17) This appliance contains batteries that are non-replaceable.

18) This appliance contains batteries that are only replaceable by skilled persons.

19) For the purposes of recharging the battery, only use the detachable supply unit provided with this appliance.

* If you feel abnormal during use, please stop using immediately and consult your doctor. Otherwise may it be harmful to health.

1. People who are susceptible to electromagnetic malfunctions such as cardiac pacemakers and implantable medical electronic devices.
2. People who are supposed to lose their sense of temperature.
3. People with malignant tumors.
4. People with heart disease.
5. A person who is restless early in pregnancy or shortly after giving birth.
6. One who has an abnormal or curved spine.
7. A person with a disturbance of perception caused by a high degree of peripheral circulation in diabetes.
8. Someone with a scar on their skin.
9. People who need quiet and people who are in poor health.
10. People with osteoporosis, severe fractures, sprain, and muscle pain, etc., acute (head) pain.

The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using this product.

Only use indoors!

Do not use the device in damp rooms (e.g. when bathing or showering)

• Do not wash!

• Do not dry clean!

Correct Disposal of this product

This marking indicates that this product should not be disposed with other household wastes throughout the EU. To prevent possible harm to the environment or human health from uncontrolled waste disposal, recycle it responsibly to promote the sustainable reuse of material resources. To return your used device, please use the return and collection systems or contact the retailer where the product was purchased. They can take this product for environmental safe recycling.

Technological parameter

Product: Handheld Massager

Model: FeelGood2 (Series JA)

Input: DC9V === 1A

Power: 9W

Massage cycle: 15 minutes

Operating time: about 50 minutes

Chargin time: about 3.5 hours

Manufacturer: XIAMEN EMOKA HEALTHSCIENCE &
TECHNOLOGY CO. LTD.

Production date: on the massager

Made in China

Manufacturer's address: D Building No. 100, JINFU ROAD,
Chengnan Industrial Zone, TONG'AN District, Xiamen, China

Phone: 0086-755-25122982

General importer in Russia:

OOO "Max Concept", OGRN 1147453008532

Moscow, Lyubertsy district, Kotelniki, 1st Pokrovsky

PR-d, 5

phone 8-800-700-14-29

www.ergonova.ru

info@ergonova.ru

Full list of representative offices: <http://ergonova.ru/about>

/stores/ all copyrights to texts and images

owned by Ergonova GmbH and may not be used otherwise

only with the consent of the copyright holder.

Copyright © 2012.

Ergonova GmbH

All rights reserved.